



Sanity 2.0 System

Fully automated and integrated molecular diagnostic platform

Features



<5 min 
of hands-on time



Maximizing Multiplexing Capacity

Empowered by patented MMCA® technology, Sanity 2.0 System is able to detect multiple targets in each fluorescent channel.



Built-in touchscreen and printer

Results can be printed out directly once complete the run.



Sample in, result out

Simple workflow effectively reduces processing time.



UV decontamination

UV lamp with a set timer to minimize the risk of contamination.



Pre-loaded reagents for easy handling

The unique design of reagent cartridges is suitable for fast-paced clinical labs.



Small footprint


Just 26 cm of bench space is needed.

Simple Workflow

1 Scan reagent barcode



2 Load sample

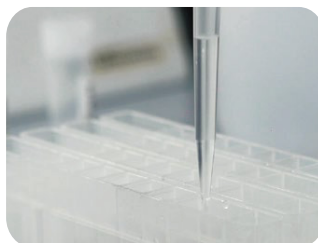

3 Insert cartridges and start the run



Tests

Tuberculosis (TB)

Cat. No.	Product Name	Description	Size
608115 608815	MTB DNA Extraction Kit (Sanity 2.0)	DNA extraction from sputum, TB culture, bronchial washing fluid, bronchoalveolar lavage fluid (BALF), etc.	24 Tests/Kit
801176*	MTB Test Kit (Sanity 2.0)	Mycobacterium Tuberculosis Complex (MTBC)	24 Tests/Kit
801177*	MDR-TB Test Kit (Sanity 2.0)	Rifampicin and isoniazid resistance mutations	24 Tests/Kit
801181*	MTB/MDR Test Kit (Sanity 2.0)	· MTBC · Rifampicin-resistance mutations · Isoniazid-resistance mutations	24 Tests/Kit
801180**	Mycobacteria Identification Kit (Sanity 2.0)	19 clinically relevant mycobacteria, including MTBC and nontuberculous mycobacteria (NTM)	24 Tests/Kit

*Use with Cat.No.608115 for extraction.

**Use with Cat.No.608815 for extraction.

Respiratory Pathogen Panel (RPP)

Cat. No.	Product Name	Description	Size
608101 608501	Sanity-2.0-Pathogen Nucleic Acid Extraction Kit	Nucleic acid extraction from saliva, BALF, sputum, swabs, etc.	24 Tests/Kit
801406**	SARS-CoV-2 Test Kit (Sanity 2.0)	Detection of SARS-CoV-2	24 Tests/Kit
801112*	Sanity-2 Respiratory Pathogen Panel (RPP4)	· Influenza A virus · RSV · Influenza B virus · SARS-CoV-2	24 Tests/Kit
801918**	Sanity-2 Respiratory Pathogen Panel (RPP40)	Detection of 40 pathogens of respiratory tract: 19 viruses · Adenovirus · Influenza B virus · Bocavirus · MERS-CoV · Coronavirus 229E/OC43/NL63/HKU1 · Parainfluenza Virus I/II/III/IV · Enterovirus · RSV A/B · Human Metapneumovirus · Rhinovirus · Influenza A virus · SARS-CoV-2 21 bacteria · <i>A. baumannii</i> · <i>C. trachomatis</i> · <i>M. catarrhals</i> · <i>S. epidermidis</i> · <i>A. fumigatus</i> · <i>Cryptococcus</i> · <i>M. pneumoniae</i> · <i>S. pneumoniae</i> · <i>B. paraptussis</i> · <i>E. coil</i> · <i>P. jirovecii</i> · <i>S. pyogenes</i> · <i>B. pertussis</i> · <i>H. influenzae</i> · <i>P. aeruginosa</i> · <i>C. pneumoniae</i> · <i>K. pneumoniae</i> · <i>Rickettsia</i> · <i>C. psittaci</i> · <i>L. pneumophila</i> · <i>S. aureus</i>	24 Tests/Kit

*Use with Cat.No.608101 for extraction. Must process 4 samples in a batch.

**Use with Cat.No.608501 for extraction. Allows for flexible loading of 1-4 samples.

HPV

Cat. No.	Product Name	Description	Size
608112	HPV DNA Extraction Kit (Sanity 2.0)	DNA extraction from cervical secretions	24 Tests/Kit
801170	High Risk HPV Genotyping Kit (Sanity 2.0)	Identification of 14 high-risk HPV types: 16, 18, 31, 33, 35, 39, 45, 51, 52, 56, 58, 59, 66, 68	24 Tests/Kit
801916	Low Risk HPV Genotyping Kit (Sanity 2.0)	Identification of 14 low-risk HPV types: 6, 11, 26, 40, 42, 43, 44, 53, 54, 61, 69, 70, 72, 81	24 Tests/Kit

STI

Cat. No.	Product Name	Description	Size
801910	STI-8 Detection kit (Sanity 2.0) For research use only	Identification of 8 STI pathogens: CT, NG, TP, TV, UU, HD, HSV-1, HSV-2	24 Tests/Kit
	STI-7 Detection Kit (Sanity 2.0) In development	Identification of 7 STI pathogens: CT, NG, UU, MG, UP, MH, HSV-2	
	STI-15 Detection Kit (Sanity 2.0) In development	Identification of 15 STI pathogens: CT, NG, TP, TV, UU, HD, MG, UP, KG, HSV-1, HSV-2, HCMV, VZV, MH, HHV-8	

Scenes



Technical Specifications

Performance Parameters	Specifications
Technology	Multicolor Melting Curve Analysis (MMCA [®]) based real-time qPCR
Sample Throughput	4 samples
Sample Volume	1 μ L~2000 μ L
Nucleic Acid Extraction Method	Magnetic bead-based
Dimension	260mm (L) x 540mm (W) x 570mm (H)
Weight	30 kg
Heating Rate	$\geq 5.0^{\circ}\text{C}/\text{s}$
Cooling Rate	$\geq 3.5^{\circ}\text{C}/\text{s}$
Temperature Accuracy	$\leq 0.1^{\circ}\text{C}$
Temperature Uniformity	$\leq \pm 0.1^{\circ}\text{C}$
Detection Channels	4 channels
Decontamination	UV light lamp
Printing function	Built-in printer
Control Panel	Touchscreen or external computer

Order Information

Cat. No.	Product Name	Size
505101	Sanity 2.0 System	1 set

Xiamen Zeesan Biotech Co., Ltd.

No. 884-1 and 884-2 Lianting Road,
(Xiang An) Industrial Area, Torch High-Tech Zone,
Xiamen City, Fujian Province, 361101, P.R. China

Phone: +86-592-7615091
E-mail: info@zsandx.com
Website: www.zeesandx.com

